



Compact Blood Bank Refrigerators





RSBG1032MD

Glass Door Compact Blood Bank

Part Number	RSBG1032MD
Capacity (Litres)	110
Capacity (250ml bags)	30 (approx)
Style	Compact
Exterior Dimensions (HxWxD)	774 x 592 x 633 mm
Interior Dimensions (HxWxD)	450 x 505 x 485 mm
Weight	59.5 kg
Energy Consumption	2.7 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Glass
Door Hinging	Right
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Shelves	2
Baskets (optional extra)	Max 6 x BABB4005
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in accordance with Directive 93/42/EEC.	✓

Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991. This blood bank is fitted with safety features such as a seven day chart recorder to keep you in the knowledge that your product is stored correctly and safely.



Hygienic touch screen
controller



Digital keypad
door lock

Fitted with a hygienic, wipe clean, touch screen controller with alarm battery back up, the clear LED display shows all alarms visually as well as emitting an audible warning. Both the air and load temperatures are monitored and there is a door open and mains failure alarm. With remote alarm contacts fitted as standard, this blood bank refrigerator can easily be linked up to most alarm systems.

We also supply a range of blood baskets purpose made for shelf models. These are designed for good air circulation around the product while keeping labels visible for ease of access and organisation.

Labcold blood banks are tested for 72 hours and calibrated by our UKAS accredited temperature calibration laboratory #8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.


RSBT1056MD Compact Blood Bank



Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991. This blood bank is fitted with safety features such as a seven day chart recorder to keep you in the knowledge that your product is stored correctly and safely.



Hygienic touch screen
controller

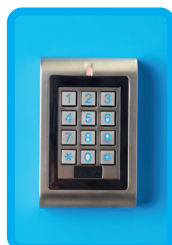

Digital keypad lock

Fitted with a hygienic, wipe clean, touch screen controller with alarm battery back up, the clear LED display shows all alarms visually as well as emitting an audible warning. Both the air and load temperatures are monitored as well as a door open and mains failure alarm. With remote alarm contacts fitted as standard, this blood bank refrigerator can be linked up to most alarm systems.

We also supply a range of blood baskets purpose made for shelf models. These are designed for good air circulation around the product while keeping labels visible for ease of access and organisation.

Labcold blood banks are tested for 72 hours and calibrated before they are dispatched and all come with a FREE 2 year parts and labour warranty.

Part Number	RSBT1056MD
Capacity (Litres)	118
Capacity (250ml bags)	56 (approx)
Style	Compact
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 545 x 450mm
Weight	77kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Solid
Door Hinging	Right
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Shelves	2
Baskets (optional extra)	Max 6
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in accordance with Directive 93/42/EEC.	✓

RSBG1056MD

Glass Door Compact Blood Bank



Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991. This blood bank is fitted with safety features such as a seven day chart recorder to keep you in the knowledge that your product is stored correctly and safely.



Hygienic touch screen controller



Digital keypad door lock

Fitted with a hygienic, wipe clean, touch screen controller with alarm battery back up, the clear LED display shows all alarms visually as well as emitting an audible warning. Both the air and load temperatures are monitored as well as a door open and mains failure alarm. With remote alarm contacts fitted as standard, this blood bank refrigerator can be linked up to most alarm systems.

We also supply a range of blood baskets purpose made for shelf models. These are designed for good air circulation around the product while keeping labels visible for ease of access and organisation.

Labcold blood banks are tested for 72 hours and calibrated by our UKAS accredited temperature laboratory #8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.

Part Number	RSBG1056MD
Capacity (Litres)	118
Capacity (250ml bags)	56 (approx)
Style	Compact
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 545 x 450mm
Weight	80kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Glass
Door Hinging	Right
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Shelves	2
Baskets (optional extra)	Max 6
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in accordance with Directive 93/42/EEC.	✓




RDBT1068MD Compact Blood Bank

Part Number	RDBT1068MD
Capacity (Litres)	118
Capacity (250ml bags)	68 (approx)
Style	Compact
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 560 x 475mm
Weight	84kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Solid
Door Hinging	Right
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Drawers	2
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓

Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991, this blood bank is fitted with safety features such as a seven day chart recorder to keep you in the knowledge that your product is stored correctly and safely.



Hygienic touch screen controller

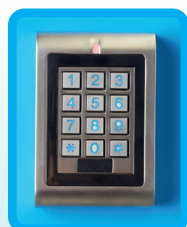

Secure digital keypad access

Fitted with a hygienic, wipe clean, touch screen controller with alarm battery back up, the clear LED display shows all alarms visually as well as emitting an audible warning. Both the air and load temperatures are monitored as well as a door open and mains failure alarm. With remote alarm contacts fitted as standard, this blood bank refrigerator can be linked up to most alarm systems.

Not only do Labcold blood banks ensure your products are stored correctly, but the digital electronic door lock provides easy access and added security. You can have all the security and safety without the need for keys.

Labcold blood banks are tested for 72 hours and calibrated in our UKAS calibration laboratory No:8898 before they are dispatched and all come with a FREE 2 year warranty.

RDBG1068MD


Compact Glass Door Blood Bank



Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991, this blood bank is fitted with safety features such as a seven day chart recorder to keep you in the knowledge that your product is stored correctly and safely.



Hygienic touch screen controller

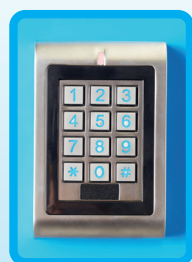

Digital keypad door lock

Fitted with a hygienic, wipe clean, touch screen controller with alarm battery back up, the clear LED display shows all alarms visually as well as emitting an audible warning. Both the air and load temperatures are monitored as well as a door open and mains failure alarm. With remote alarm contacts fitted as standard, this blood bank refrigerator can be linked up to most alarm systems.

Labcold blood banks are tested for 72 hours and calibrated by our UKAS accredited temperature calibration laboratory #8898 before they are dispatched and all come with a FREE 2 year warranty.

Part Number	RDBG1068MD
Capacity (Litres)	118
Capacity (250ml bags)	68 (approx)
Style	Compact
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 560 x 475mm
Weight	85kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Glass
Door Hinging	Right
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Drawers	2
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓




RSBT2056MD Compact Dual Compressor Blood Bank

Part Number	RSBT2056MD
Capacity (Litres)	118
Capacity (250ml bags)	56 (approx)
Style	Freestanding
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 545 x 450mm
Weight	77kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Solid
Door Hinging	Right
Dual Refrigeration	✓
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Shelves	2
Baskets (optional extra)	Max 6
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓

Ideal for A&E, Maternity, Community Hospitals, and Hospices etc, fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.



Hygienic touch screen controller

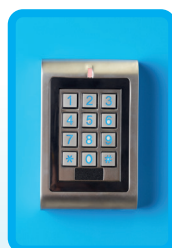

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified.

We also supply a range of blood baskets purpose made for shelf models. These are designed for good air circulation around the product while keeping labels visible for ease of access and organisation.

Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991.

Labcold blood banks are tested for 72 hours and calibrated by our UKAS temperature calibration laboratory #8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.




LABCOLD™
medical & scientific refrigeration



Ideal for A&E, Maternity, Community Hospitals, and Hospices etc, fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.



Hygienic touch screen controller

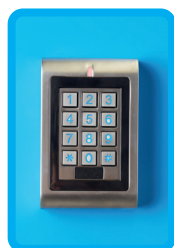

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified. We also supply a range of blood baskets purpose made for shelf models. These are designed for good air circulation around the product while keeping labels visible for ease of access and organisation. Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991.

Labcold blood banks are tested for 72 hours and calibrated before they are dispatched and all come with a FREE 2 year parts and labour warranty.

RSBG2056MD Compact Dual Compressor Blood Bank

Part Number	RSBG2056MD
Capacity (Litres)	118
Capacity (250ml bags)	56 (approx)
Style	Freestanding
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 545 x 450mm
Weight	80kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Glass
Door Hinging	Right
Dual Refrigeration	✓
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back Up	✓
Shelves	2
Baskets (optional extra)	Max 6
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓

01256 705 570

www.labcold.com



RDBT2068MD


Compact Dual Compressor Blood Bank



Ideal for A&E, Maternity, Community Hospitals, and Hospices etc, fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.



Hygienic touch screen controller

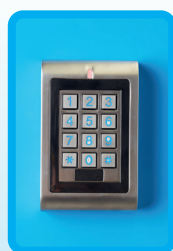

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified. Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991.


Labcold blood banks are tested for 72 hours and calibrated before they are dispatched and all come with a FREE 2 year parts and labour warranty.

Part Number	RDBT2068MD
Capacity (Litres)	118
Capacity (250ml bags)	68 (approx)
Style	Freestanding
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 560 x 475mm
Weight	84kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Solid
Door Hinging	Right
Dual Refrigeration	✓
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back up (Display)	✓
Drawers	2
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓


RDBG2068MD Compact Dual Compressor Blood Bank



Ideal for A&E, Maternity, Community Hospitals, and Hospices etc, fitted with a Dual Compressor to ensure greater temperature control and peace of mind in the event of failure of one of the systems. Each complete refrigeration system provides 100% performance, sufficient to cool down and maintain temperatures correctly without help from the second compressor.

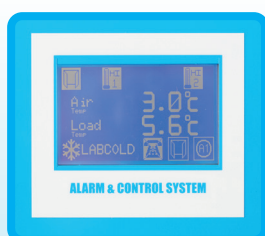

Hygienic touch screen
alarm controller

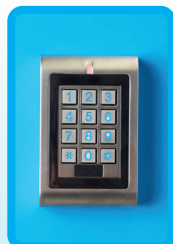

Digital keypad lock

An asymmetric change over ensures even distribution of workload, resulting in even wear of the compressors and also allows for one to defrost while the other is still running. The two independent systems provide back up should an issue develop plus there is no need to move contents immediately if a problem is identified. Designed and manufactured in the UK in our Basingstoke factory, this Blood Bank refrigerator is a certified Medical Device in accordance with Directive 93/42/EEC and is compliant with BS4376:Part 1:1991.

Labcold blood banks are tested for 72 hours and calibrated in our UKAS calibration laboratory No:8898 before they are dispatched and all come with a FREE 2 year parts and labour warranty.

Part Number	RDBG2068MD
Capacity (Litres)	118
Capacity (250ml bags)	68 (approx)
Style	Freestanding
Exterior Dimensions (HxWxD)	1040 x 645 x 680mm
Interior Dimensions (HxWxD)	445 x 560 x 475mm
Weight	85kg
Energy Consumption	2.4 kWh/24h
Operating Temperature	+4°C +/- 2°C
Interior Construction	Stainless Steel
Exterior Construction	Powder Coated Steel
Door Style	Glass
Door Hinging	Right
Dual Refrigeration	✓
Open Door alarm	✓
Load Temperature Alarm	✓
Air Temperature Alarm	✓
Mains Failure Alarm	✓
Air Temperature Display	✓
Load Temperature Display	✓
7 Day Chart Recorder	✓
Lock	Digital Keypad
Battery Back up (Display)	✓
Drawers	2
Remote alarm contacts	✓
Access Port	✓
Auto Defrost	✓
Internal Light	✓
Eco Friendly Refrigerant	✓
Certified Medical Device in Accordance with Directive 93/42/EEC.	✓